

POLYCLONAL ANTIBODY

Anti-HA-tag

Code No.
561-5

Quantity
500 µL

Form
Rabbit IgG

BACKGROUND: Epitope tagging has widely been accepted technique that fuse an epitope peptide to a certain protein as a marker for gene expression. With this technique, the gene expression can be easily monitored on western blotting, immunoprecipitation and immunofluorescence utilizing with an antibody that recognizes such an epitope. Amino acid sequences that are widely used for the epitope tagging are as follow; YPYDVPDYA (HA-tag), EQKLISEEDL (Myc-tag) and YTDIEMNRLGK (VSV-G-tag), which corresponding to the partial peptide of Influenza hemagglutinin protein, Human c-myc gene product and Vesicular stomatitis virus glycoprotein respectively.

SOURCE: This antibody was purified from rabbit serum using protein A agarose. The rabbit was immunized with KLH conjugated synthetic peptide, KLH-YPYDVPDYA (HA-tag).

FORMULATION: 500 µL volume of PBS containing 50% glycerol, pH 7.2. No preservative is contained.

STORAGE: This antibody solution is stable for one year from the date of purchase when stored at -20°C.

REACTIVITY: This antibody reacts with HA-tag on Western blotting.

APPLICATIONS:

Western blotting: 1:1,000 for chemiluminescence detection system

Immunoprecipitation: 1 µL

Immunohistochemistry: Not tested

Immunocytochemistry: 1:200

Flow cytometry: Not tested

Detailed procedure is provided in the following **PROTOCOLS**.

INTENDED USE:


For Research Use Only. Not for use in diagnostic procedures.

REFERENCES:

- 1) Fukunaka, A., *et al.*, *J. Biol. Chem.* **286**, 16363-16373 (2011)
- 2) Mimura, Y., *et al.*, *J. Cell Sci.* **123**, 2014-2024 (2010)
- 3) Murano, K., *et al.*, *Mol. Cell Biol.* **28**, 3114-3126 (2008)
- 4) Shirai, A., *et al.*, *J. Biol. Chem.* **283**, 10745-10752 (2008)

- 5) Sun, P., *et al.*, *Genes Cells* **11**, 1097-1113 (2006)
- 6) Arao, Y., *et al.*, *Mol. Endocrinol* **18**, 2255-2267 (2004)
- 7) Myojin, R., *et al.*, *Biol. Reprod* **71**, 926-932 (2004)
- 8) Zhang, A., *et al.*, *J. Biol. Chem.* **279**, 33799-33805 (2004)
- 9) Murayama, Y., *et al.*, *J. Cell Sci.* **117**, 3379-3388 (2004)
- 10) Saito, T. T., *et al.*, *Nucleic Acids Res.* **32**, 3325-3339 (2004)
- 11) Toji, S., *et al.*, *Genes Cells* **9**, 383-397 (2004)
- 12) Yasuda, S., *et al.*, *Endocrinology* **145**, 1730-1738 (2004)
- 13) Sekiguchi, T., *et al.*, *J. Biol. Chem.* **279**, 8343-8350 (2004)
- 14) Abe, T., *et al.*, *J. Cell Sci.* **116**, 155-168 (2003)
- 15) Sun, P., *et al.*, *J. Biol. Chem.* **278**, 4063-4071 (2003)
- 16) Furusawa, M., *et al.*, *J. Biol. Chem.* **277**, 50885-50892 (2002)
- 17) Katsu, R., *et al.*, *J. Biol. Chem.* **277**, 44220-44228 (2002)
- 18) Zama, T., *et al.*, *J. Biol. Chem.* **277**, 23909-23918 (2002)
- 19) Nishizawa, H., *et al.*, *J. Biol. Chem.* **277**, 1586-1592 (2002)
- 20) Kamimoto, T., *et al.*, *J. Biol. Chem.* **276**, 37520-37528 (2001)
- 21) Kojima, T., *et al.*, *J. Biol. Chem.* **276**, 32247-32256 (2001)
- 22) Sekiguchi, T., *et al.*, *J. Biol. Chem.* **276**, 7246-7257 (2001)
- 23) Izumi, M., *et al.*, *Nucleic Acid Res.* **28**, 4769-4777 (2000)
- 24) Matsuda, S., *et al.*, *J. Biol. Chem.* **274**, 34515-34518 (1999)
- 25) Takahashi, M., *et al.*, *J. Virol.* **73**, 7812-7816 (1999)
- 26) Koizumi, J., *et al.*, *J. Biol. Chem.* **274**, 11125-11131 (1999)
- 27) Kondo, T., *et al.*, *Mol. Cell Biol.* **19**, 1136-1143 (1999)
- 28) Furuse, M., *et al.*, *J. Cell Biol.* **143**, 391-401 (1998)
- 29) Shimomura, T., *et al.*, *Mol. Cell Biol.* **18**, 5485-5491 (1998)

This antibody is used in these references.

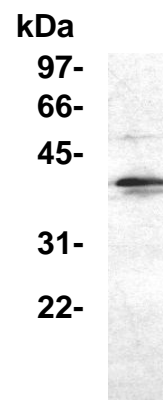


Western blot analysis of HA-tagged I-κB protein expression in 293T using 561-5.

PROTOCOLS:

SDS-PAGE & Western Blotting

- 1) Mix the sample with equal volume of Laemmli's sample buffer.
- 2) Boil the samples for 2 minutes and centrifuge. Load 10 μ L of the sample per lane in a 1-mm-thick SDS-polyacrylamide gel for electrophoresis.
- 3) Blot the protein to a polyvinylidene difluoride (PVDF) membrane at 1 mA/cm² for 1 hour in a semi-dry transfer system (Transfer Buffer: 25 mM Tris, 190 mM glycine, 20% MeOH). See the manufacturer's manual for precise transfer procedure.
- 4) To reduce nonspecific binding, soak the membrane in 10% skimmed milk (in PBS, pH 7.2) for 1 hour at room temperature, or overnight at 4 °C.
- 5) Incubate the membrane with primary antibody diluted with PBS, pH 7.2 containing 1% skimmed milk as suggest in the **APPLICATIONS** for 1 hour at room temperature. (The concentration of antibody to be used will be depend on condition.)
- 6) Wash the membrane with PBS-T [0.05% Tween-20 in PBS] (5 minutes x 3 times).
- 7) Incubate the membrane with the 1:10,000 HRP-conjugated anti-rabbit IgG (MBL; code no. 458) diluted with 1% skimmed milk (in PBS, pH 7.2) for 1 hour at room temperature.
- 8) Wash the membrane with PBS-T (5 minutes x 6 times).
- 9) Wipe excess buffer on the membrane, then incubate it with appropriate chemiluminescence reagent for 1 minute. Remove extra reagent from the membrane by dabbing with paper towel, and seal it in plastic wrap.
- 10) Expose to an X-ray film in a dark room for 5 minutes. Develop the film as usual. The condition for exposure and development may vary.

kDa

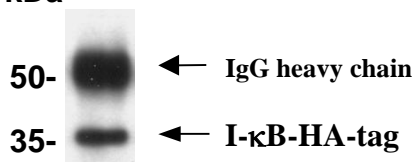

Immunoprecipitation of HA-tagged I-κB protein from transfectant with 561-5. After immunoprecipitated with the antibody, immunocomplex was resolved on SDS-PAGE and immunoblotted with 561-5.

Immunoprecipitation

- 1) Wash the cells 3 times with PBS and suspend with 10 volume of cold Lysis buffer (50 mM Tris-HCl pH 7.2, 250 mM NaCl, 0.1% NP-40, 2 mM EDTA, 10% glycerol) containing appropriate protease inhibitors. Incubate it at 4 °C with rotating for 30 minutes, then sonicate briefly (up to 10 seconds).
- 2) Centrifuge the tube at 12,000 x g for 10 minutes at 4 °C and transfer the supernatant to another tube.
- 3) Add primary antibody as suggest in the

APPLICATIONS into 200 μ L of the supernatant. Mix well and incubate with gentle agitation for 30-120 minutes at 4 °C. Add 20 μ L of 50% protein A agarose beads resuspended in the cold Lysis buffer. Mix well and incubate with gentle agitation for 60 minutes at 4 °C.

- 4) Wash the beads 3-5 times with the cold Lysis buffer (centrifuge the tube at 2,500 x g for 10 seconds).
- 5) Resuspend the beads in 20 μ L of Laemmli's sample buffer, boil for 3-5 minutes, and centrifuge for 5 minutes. Use 10 μ L/lane for the SDS-PAGE analysis. (See **SDS-PAGE & Western blotting.**)

Immunocytochemistry

- 1) Culture the cells in the appropriate condition on a glass slide.
- 2) Fix the cells by immersing the slide in acetone for 10 minutes on ice.
- 3) Air dry the slides.
- 4) Add the primary antibody diluted with PBS as suggest in the **APPLICATIONS** onto the cells and incubate for 30 minutes at room temperature. (Optimization of antibody concentration or incubation condition are recommended if necessary.)
- 5) Prepare a wash container such as a 500 mL beaker with a stirrer. Then wash the cells on the glass slide by soaking the slide with a plenty of PBS in the wash container for 5 minutes. Take care not to touch the cells. Repeat another wash once more.
- 6) Add 30 μ L of 1:40 FITC conjugated anti-rabbit IgG (MBL; code no. IM-0833) diluted with PBS onto the cells. Incubate for 30 minutes at room temperature. Keep out light by aluminum foil.
- 7) Wash the slide in a plenty of PBS as in the step 5).
- 8) Wipe excess liquid from slide but take care not to touch the cells. Never leave the cells to dry.
- 9) Promptly add mounting medium onto the slide, then put a cover slip on it.

RELATED PRODUCTS

Antibodies

561	anti-HA-tag (polyclonal) (0.1 mL)
561-5	anti-HA-tag (polyclonal) (0.5 mL)
561-7	anti-HA-tag HRP-Direct (polyclonal)
561-8	anti-HA-tag-agarose (polyclonal)
M180-3	anti-HA-tag (TANA2) (200 μ L)
M180-3S	anti-HA-tag (TANA2) (50 μ L)
M180-A48	anti-HA-tag-Alexa Fluor [®] 488 (TANA2)
M180-A59	anti-HA-tag-Alexa Fluor [®] 594 (TANA2)
M180-A64	anti-HA-tag-Alexa Fluor [®] 647 (TANA2)
M132-3	anti-HA-tag (5D8)
M185-3L	anti-DDDDK-tag (FLA-1) (1 mL)
M185-3LL	anti-DDDDK-tag (FLA-1) (5 mL)
M185-3S	anti-DDDDK-tag (FLA-1) (50 μ L)
M185-A48	anti-DDDDK-tag-Alexa Fluor [®] 488 (FLA-1)
M185-A59	anti-DDDDK-tag-Alexa Fluor [®] 594 (FLA-1)
M185-A64	anti-DDDDK-tag-Alexa Fluor [®] 647 (FLA-1)
PM020	anti-DDDDK-tag (polyclonal)
PM020-7	anti-DDDDK-tag HRP-Direct (polyclonal)

PM020-8 anti-DDDDK-tag-agarose (polyclonal)
D291-3 anti-His-tag (OGHis) (200 μ L)
D291-3S anti-His-tag (OGHis) (50 μ L)
D291-6 anti-His-tag-biotin (OGHis)
D291-7 anti-His-tag HRP-Direct (OGHis)
D291-8 anti-His-tag-agarose (OGHis)
D291-A48 anti-His-tag-Alexa Fluor[®] 488 (OGHis)
D291-A59 anti-His-tag-Alexa Fluor[®] 594 (OGHis)
D291-A64 anti-His-tag-Alexa Fluor[®] 647 (OGHis)
M089-3 anti-His-tag (6C4)
M136-3 anti-His-tag (2D8)
PM032 anti-His-tag (polyclonal)
PM032-8 anti-His-tag-agarose (polyclonal)
598 anti-GFP (polyclonal)
598-7 anti-GFP HRP-Direct (polyclonal)
M048-3 anti-GFP (1E4)
D153-3 anti-GFP (RQ2)
D153-A48 anti-GFP-Alexa Fluor[®] 488 (RQ2)
D153-A59 anti-GFP-Alexa Fluor[®] 594 (RQ2)
D153-A64 anti-GFP-Alexa Fluor[®] 647 (RQ2)
D153-8 anti-GFP-agarose (RQ2)
PM005 anti-RFP (polyclonal)
PM005-7 anti-RFP HRP-Direct (polyclonal)
M155-3 anti-RFP (8D6)
M165-3 anti-RFP (3G5)
M165-8 anti-RFP-agarose (3G5)
M071-3 anti-GST-tag (3B2)
PM013 anti-GST-tag (polyclonal)
PM013-7 anti-GST-tag HRP-Direct (polyclonal)
M013-3 anti-Thioredoxin (2C9)
M094-3 anti- β -galactosidase (5A3)
PM049 anti- β -galactosidase (polyclonal)
M095-3 anti-Luciferase (2D4)
PM016 anti-Luciferase (polyclonal)
PM047 anti-Renilla Luciferase (polyclonal)
M091-3 anti-MBP (1G12)
PM015 anti-Chitin Binding Domain (polyclonal)
M047-3 anti-Myc-tag (PL14)
M047-6 anti-Myc-tag-biotin (PL14)
M047-7 anti-Myc-tag HRP-Direct (PL14)
M047-8 anti-Myc-tag-agarose (PL14)
M047-A48 anti-Myc-tag-Alexa Fluor[®] 488 (PL14)
M047-A59 anti-Myc-tag-Alexa Fluor[®] 594 (PL14)
M047-A64 anti-Myc-tag-Alexa Fluor[®] 647 (PL14)
562 anti-Myc-tag (polyclonal) (0.1 mL)
562-5 anti-Myc-tag (polyclonal) (0.5 mL)
PM021 anti-S-tag (polyclonal)
PM021-8 anti-S-tag-agarose (polyclonal)
PM022 anti-T7-tag (polyclonal)
PM022-8 anti-T7-tag-agarose (polyclonal)
M167-3 anti-V5-tag (1H6)
PM003 anti-V5-tag (polyclonal)
PM003-7 anti-V5-tag HRP-Direct (polyclonal)
PM003-8 anti-V5-tag-agarose (polyclonal)
563 anti-VSV-G-tag (polyclonal)
563-8 anti-VSV-G-tag-agarose (polyclonal)

Smart-IP series

3190 Magnetic Rack
M132-9 anti-HA-tag-magnetic beads (5D8)
M185-9 anti-DDDDK-tag-magnetic beads (FLA-1)
D291-3 anti-His-tag-magnetic beads (OGHis)
D153-9 anti-GFP-magnetic beads (RQ2)
M165-9 anti-RFP-magnetic beads (3G5)
M047-9 anti-Myc-tag-magnetic beads (PL14)
M167-9 anti-V5-tag-magnetic beads (1H6)

Protein Purification Kit

3305 c-Myc-tagged Protein MILD PURIFICATION KIT
3305A c-Myc-tagged Protein MILD PURIFICATION KIT (Trial Kit)
3306 c-Myc-tagged Protein MILD PURIFICATION GEL (1 mL gel, 1 mg peptide)
3307 c-Myc-tagged Protein MILD PURIFICATION GEL (5 mL gel, 5 mg peptide)
3300-205 c-Myc tag peptide (5 mg)
3310 His-tagged Protein PURIFICATION KIT
3310A His-tagged Protein PURIFICATION KIT (Trial Kit)
3310-205 His-tag peptide (10mg)
3311 His-tagged Protein PURIFICATION GEL (1 mL gel, 10 mg peptide)
3312 His-tagged Protein PURIFICATION GEL (5 mL gel, 50 mg peptide)
3315 V5-tagged Protein PURIFICATION KIT
3315A V5-tagged Protein PURIFICATION KIT (Trial Kit)
3320 HA-tagged Protein PURIFICATION KIT
3320A HA-tagged Protein PURIFICATION KIT (Trial Kit)
3325 DDDDK-tagged Protein PURIFICATION KIT
3325A DDDDK-tagged Protein PURIFICATION KIT (Trial Kit)
3325-205 DDDDK-tag peptide (5 mg)
3326 DDDDK-tagged Protein PURIFICATION GEL with Elution Peptide (1 mL gel, 5 mg peptide)
3327 DDDDK-tagged Protein PURIFICATION GEL with Elution Peptide (5 mL gel, 25 mg peptide)
3328 DDDDK-tagged Protein PURIFICATION GEL (5 mL gel)
3329 DDDDK-tagged Protein PURIFICATION GEL (25 mL gel)

Other related antibodies and kits are also available.
Please visit our website at <https://ruo.mbl.co.jp/>